

Aspect100EP


PRODUCT DATA

Features:

- 61" combined unobstructed tabletop viewing
- Motorized angulation: 75/15°
- Adjustable height: 29" to 39"
- Table dimensions: 22" x 84"
- Capacity: 350#
- Table weight: ~475#
- Low attenuation carbon fiber tabletop
- Motorized elevation with hand control
- Locking swivel casters
- 1/2" table pad and restraint straps
- Footboard
- Hand holds

Accessories:

- Tableside radiation shield
- Footswitch

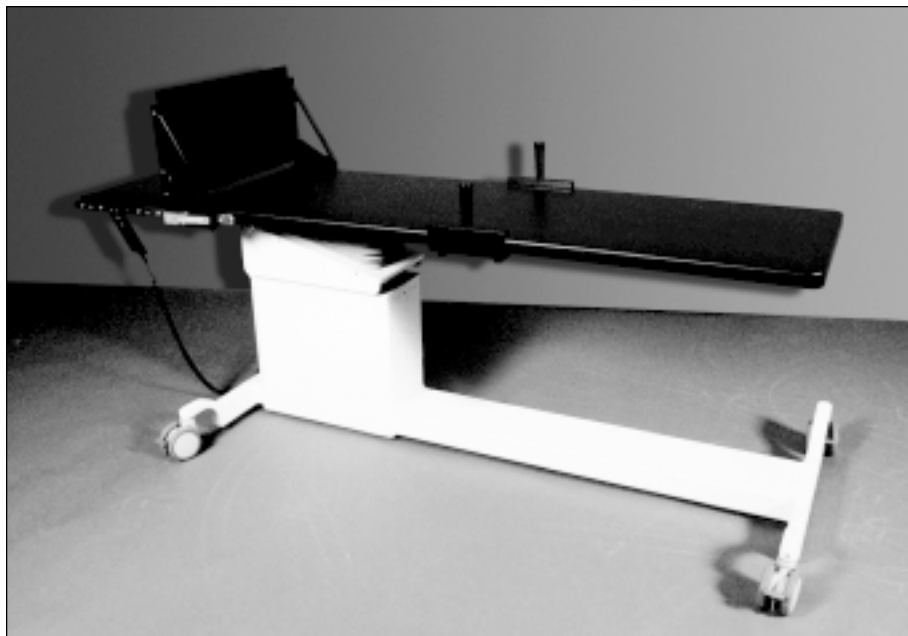

Aspect100EP Fluoroscopic Examining Table

Image Diagnostics' *Aspect100EP* is specifically designed for ease of use and high quality image outcome. *Aspect* tables benefit from a trim, lightweight design, ensuring a single operator can position the table with ease. Features unique to the *Aspect* table series provide the ultimate in table maneuverability, operator access, and patient comfort. The Image Diagnostics *Aspect100EP* is designed to meet the imaging requirements of electrophysiologic examinations.

The innovative, cantilevered table design of the *Aspect100EP* significantly enhances mobile C-arm positioning and provides visualization capabilities free from obstructions or artifacts. A unique, center-rib base design allows the operator to stand tableside without having to step on or around casters. Approaches to the table can be made from either side, as well as the head-end, allowing enhanced patient access with little or no need for repositioning.

The *Aspect100EP* is designed with an extended base that eliminates the need for heavy ballast for stability and helps protect the tabletop from impacts. Vertical travel capabilities of the *Aspect100EP* allow patient placement at a comfortable working height and facilitate its successful integration into a variety of mobile C-arm and fixed-base imaging products.

A variety of table accessories are available for patient positioning and operator radiation protection.

Aspect100EP Fluoroscopic Examining Table

FEATURES AND BENEFITS


- Carbon fiber tabletop is 84" in length with 61" imaging area; for most patients, provides total body coverage with minimal need for C-arm repositioning.
- Table width is 22"; provides excellent access for lateral C-arm positioning.
- Power-driven table angulation covers 75/15°; enhances examination capabilities and improves test accuracy.
- Motorized elevation with height range of 39" to 29" is one of the lowest available retracted heights; facilitates patient loading and unloading.

For a complete listing of Image Diagnostics' tables and products, visit our website at www.imagediagnostics.com.

Control Accessories



Footswitch

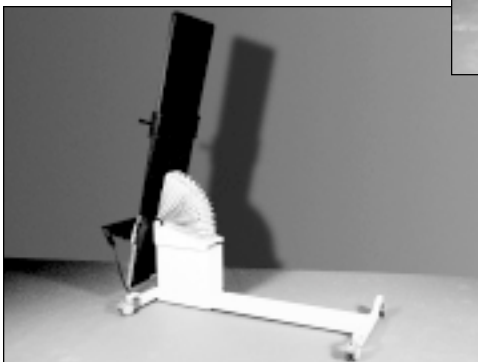

*All Aspect series tables and accessories are backed by a three year warranty on parts.
Image Diagnostics, Inc. reserves the right to update product specifications, without notification, as required.
www.imagediagnostics.com*

100EP